

# lifetime limited warranty

Your Drive branded product is warrantied to be free of defects in materials and workmanship for the lifetime of the product for the original consumer purchaser.

This device was built to exacting standards and carefully inspected prior to shipment. This Lifetime Limited Warranty is an expression of our confidence in the materials and workmanship of our products and our assurance to the consumer of years of dependable service.

This warranty does not cover device failure due to owner misuse or negligence, or normal wear and tear. The warranty does not extend to non-durable components, such as rubber accessories, casters, and grips, which are subject to normal wear and need periodic replacement.

If you have a question about your Drive device or this warranty, please contact an authorized Drive dealer.



© 2015 Medical Depot, Inc. All rights reserved.

Drive is a trademark of Medical Depot, Inc.

Port Washington N.Y. 11050 USA Made in China

Rev.5.12.02.15

exercise peddler  
with handle

item #10274



EU Authorized Representative  
Drive Medical LTD  
Ainley's Industrial Estate  
Elland, West Yorkshire,  
United Kingdom HX5 9JP


[www.drivemedical.com](http://www.drivemedical.com)

drive<sup>TM</sup>

## assembly instructions

1. Position handle with the peddler base so mounting screw holes align. Insert bolts and secure in place with wing nuts, turning wing nuts clockwise until tight.
2. **Handle Height Adjustment:** Extend handle to desired height. Insert locking collar through hole to lock handle into place.

## resistance adjustment

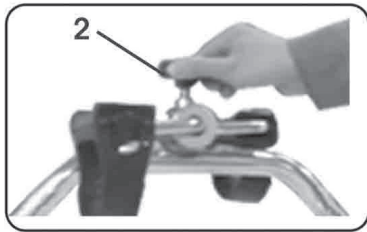

Turn the tension knob to adjust the paddle resistance. “High” means heavy resistance, “Low” means light resistance.

**WARNING**—DEVICE MAY BECOME HOT WITH EXTENDED USE.

