



#1 in Knee Bracing

CCMI Mark III

DJO Global
Vista Manufacturing Center
3151 Scott Street
Vista, CA 92081-8553
USA

Tel: (800) 336-6569
www.djoglobal.com

Fax To: (800) 457-4221

Email To: defiance.help@djoglobal.com

BILL TO:

Account Name: _____

Account #: _____ P.O. #: _____

Contact: _____ Phone: _____ Ext. _____

Address: _____

City: _____ State: _____ Zip: _____

SHIPPING INFORMATION: Same as "BILL TO"

Name: _____

Address: _____

City: _____ State: _____ Zip: _____

Next Day Air 2nd Day Delivery

3rd Day Delivery Ground Other _____

PRODUCT SELECTION: CCMI Mark III

Frame: DEFIANCE® Football/Collision

Hinge:

Elbow (CI only)

Strapping Sequence:

CI

Measurement taken by:

Cast

CCMI Mark III

Recommended Football Options

- Chamois Liners/Pads
- Reinforced Carbon Composite Frame
- Metal Caps — Power Caps
- Square Strap Ends
- Player Last Name, First Initial

- School/Team Color: _____

FOR DJO USE ONLY:

Order #: _____

Brace #: _____

PATIENT INFORMATION:


Name: _____

Age: _____ Height: _____ Weight: _____

Measurement: Right Left Reverse

Measured By: _____

New Brace Remeasurement/Repair Refurbish

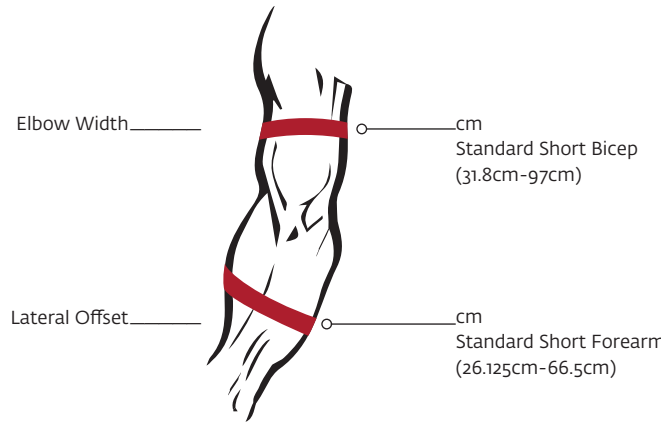

NON-METALLIC COLORS:

- Matte Black Neon Pink Neon Yellow Neon Orange
- Black Red Purple Pink
- White Blue Light Blue Yellow
- Green Orange Tan Red, White, Blue

METALLIC COLORS:

- Black Red Purple Pink
- Navy Dark Blue Bright Blue Gold
- Green Orange White Silver

MULTIPLE COLORS:

Thigh Frame: _____

Calf Frame: _____

Fade (non-metallic colors only): Yes No

GRAPHICS (extra charge):

- Tiger Red Tiger Orange Carbon Fiber
- Snake Camouflage
- Stars & Stripes (non-metallic) Stars & Stripes (metallic)

FEATURES:

- Suede Liners/Pads
- Chamois Liners/Pads
- Swiveling Strap Tabs*
- Metal Caps - PowerCaps (recommended for contact sports)
- Reinforced Carbon Composite Frame
- Installed Extension Stop: 0° 10° 20° 30° 40° (circle one)
- Installed Flexion Stop: 45° 60° 75° 90° (circle one)
- Neoprene Suspension Strap Kit

ACCESSORIES (extra charge):

- High Strength Stainless Steel Gears Air Condyle Kit
- Full Extension Lockout Stop Floam Condyle Pads
- Lycra® Suspension Sleeve Atrophy Pad Kit
- Neoprene Suspension Sleeve
- Sports Brace Cover
- Sili-Grip Strap Pad Kit
- Condyle Pad Kit


For further instructions, please refer to the CCMI Mark III instruction video and instruction card included with the instrument.

CAUTION!

FEDERAL LAW (U.S.A.) RESTRICTS THIS DEVICE TO SALE BY THE ORDER OF A LICENSED HEALTH CARE PROFESSIONAL

Custom Contour Measuring Instrument - Mark III Use Instructions

NOTICE: For optimum long term fit of the DEFIANCE® or DEFIANCE® III, measurements should be taken when there is no edema or atrophy present.

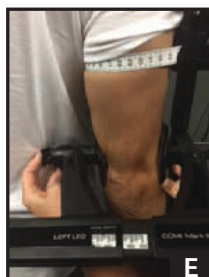

APPLICATION INFORMATION:

1. To unfold, distract and rotate swing arms of the CCMI Mark III device. **Figures A1 and A2.**
2. Have the patient stand with their back facing the practitioner, and the shoulder of the affected arm slightly extended. The Elbow should be relaxed and the palm should be facing forward. Do not allow the elbow to be flexed beyond 15 degrees. Apply device directly to bare arm. **Figure B.**
3. Turn the instrument over so that it reads opposite of the limb being measured. Pull apart the knee width caliper and position on patient's elbow. Arrows on condyle shells should align with the top of the Olecranon process and midline with the Humerous. **Figure C.**
4. For short thigh (bicep measurement) and calf (forearm measurement), extend tape measure marked in centimeters located proximal to elbow. **Figure D.**
5. To establish proper compression for elbow width measurement, support lateral side of instrument and press load plunger on medial condyle until triangles are aligned. **Figure E.**
6. Record elbow width on order form. **Figure F.**
7. Record lateral offset on order form. **Figure G.**
8. Record bicep and forearm circumference on order form. **Figure H.**